

Refrigerator/Freezer Combination LabStar Sirius

- Forced air circulation for refrigerator and freezer with automatic defrost and condensate evaporation
- 2 separate cooling circuits for refrigerator and freezer
- Cooling unit top
- 2 operating panels with:
 - Main switch
 - Electronic control with digital display
 - Visual/audible alarm device
 - Alarm warning for excess and insufficient temperature
- Door alarm
- No-volt remote alarm contact for connection to a remote alarm system
- RS 485 port
- Self-closing and lockable doors
- Door hinges on right, changeable
- Door seal changeable
- 4 castors, 2 lockable
- 4 (2x2) adjustable grid inserts in tilt-resistant guide rails
- Interior made of stainless steel CNS



LSSI 7030KEEU
(optional with feet)

Model LabStar Sirius	LSSI 3030	LSSI 7030
Version	KEWU/KEEU	KEWU/KEEU

Refrigerator:

Temperature range	°C: +1/+10	+1/+10
Content	litres (net): 160	290
Int. dimensions (WxDxH)	mm: 400x560x700	600x685x700
Grid inserts	: 2	2
Refrigerant CFC- and HCFC-free	: R134a	R134a

Freezer:

Temperature range	°C: -10/-30	-10/-30
Content	litres (net): 160	290
Int. dimensions (WxDxH)	mm: 400x560x700	600x685x700
Grid inserts	: 2	2
Refrigerant CFC- and HCFC-free	: R404A	R404A

Ext. dimensions (WxDxH)	mm: 520x700x1997	720x860x1997
Depth with door opened	mm: 1182	1540
Weight	kgs: 120	151
Mains connection	: 1/N/PE ~ 230 V 50 Hz	1/N/PE ~ 230 V 50 Hz
Shipping dimensions (WxDxH)	mm: 610x840x2200	800x900x2200
Shipping weight	kgs: 130	163

LabStar Sirius + Accessories:	Outside/ Inside	Cat. No.
LabStar Sirius LSSI 3030KEWU	white / CNS	GSLSSI3030KEWU
LabStar Sirius LSSI 3030KEEU	CNS / CNS	GSLSSI3030KEEU
LabStar Sirius LSSI 7030KEWU	white / CNS	GSLSSI7030KEWU
LabStar Sirius LSSI 7030KEEU	CNS / CNS	GSLSSI7030KEEU

ProfiTec Temperature Monitoring and Recording System (see sep. brochure). GPTEC4A

- ProfiCaller telephone dialler (see sep. brochure). on request
- Pt-100 temperature sensor, add. installed for ext. temperature recording. ZXX2-020-2L
- Grid inserts, 330x545 mm (WxD). ZLSSI4-330x545
- Grid inserts, 530x655 mm (WxD). ZLSSI4-530x655
- Guide rails for grid inserts LSSI 3030.. (1 pair). ZLSSI4-023
- Guide rails for grid inserts LSSI 7030.. (1 pair). ZLSSI4-024
- Drawer with telescopic slide-out for LSSI 3030.., made of CNS. ZLSSI4-003
- Drawer with telescopic slide-out for LSSI 7030.. made of CNS. ZLSSI4-002
- Interior socket, tropical design. ZLSSI1-001-01
- Feet instead of castors, made of CNS (50 mm height adjustable). ZLSSI5-021-125
- Duct to storage compartment 26 mm in diameter. ZLSSI1-002-26
- Door hinges factory-mounted on left. ZLSSI5-001
- MIN/MAX thermometer to be placed in refrigerators or freezers. ZXX2-019-001
- Measuring range -50/+50 °C, an add. temperature sensor is not required.
- MIN/MAX digital thermometer. Hand-held unit with digital display as well as data-hold and relative measuring function. Measuring range -200/+1370 °C. ZXX2-019-011
- An add. temperature sensor is required, eg. Cat.-No. ZXX2-022-100.
- Temperature sensor for MIN/MAX thermometer ZXX2-019-011. ZXX2-022-100
- Measuring range -100/+1000 °C. Length 2 m. Other ranges on request.

Innovation prize


ProfiTec Intelligent Temperature Monitoring and Recording System



- For up to 4 refrigerators and/or freezers
- Readings are stored over a period up to 781 days
- Warning for excess and insufficient temperature
- Includes No-volt remote alarm contact
- The measuring protocol can be fed to PC via RS 232 port

Please ask for a personal fitting option.

Recommendations for additional storage sets can be found in our separate catalogue.

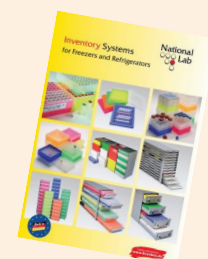

<http://www.catalopedia.de>



<http://www.national-lab.com>

National Lab GmbH
 Grambeker Weg 157
 D-23879 Mölln

Tel.: +49 45 42 / 84 91 - 70
 Fax: +49 45 42 / 84 91 - 71
 e-Mail: LabSales@NationalLab.com